TELEMAX

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Warning!

Please read the instructions carefully before using the product.

Store them in a safe place for future reference.

Safety Instructions

This device is intended for home use only.

Remove the batteries from the product if you do not intend to use it for an extended period.

Keep this device away from children and individuals with special needs.

Do not move the device while it is in use.

Do not exceed the device's weight capacity.

Do not use the device if your feet are wet.

Do not store the device in an upright position.

Do not use the device if the glass is cracked or broken.

Improper use of the product can be dangerous, making the warranty void and absolving any liability.

Telexmax disclaims any liability for damage that may occur to people, animals, or property due to non-compliance with these warnings.

Make sure that the safety instructions have been read from the start. The manual should be followed at all times during use. The scale allows for a maximum weight (load) of 180kg.

Follow the instructions below to start this device.


Product Specification

- **Size:** 300 x 260 x 18 mm
- **Power:** 2 x AAA batteries (not included)
- **Display:** LCD readout
- **Maximum Capacity:** 180 kg / 400 lb / 28 st
- **Division:** 100 g / 0.2 lb
- **Minimum Effective Weight:** 6 kg

Product Features

- Automatic on/off function
- Auto zero reset
- Weight units: kg - lb - st
- Low battery and overload indicators

Product Overview

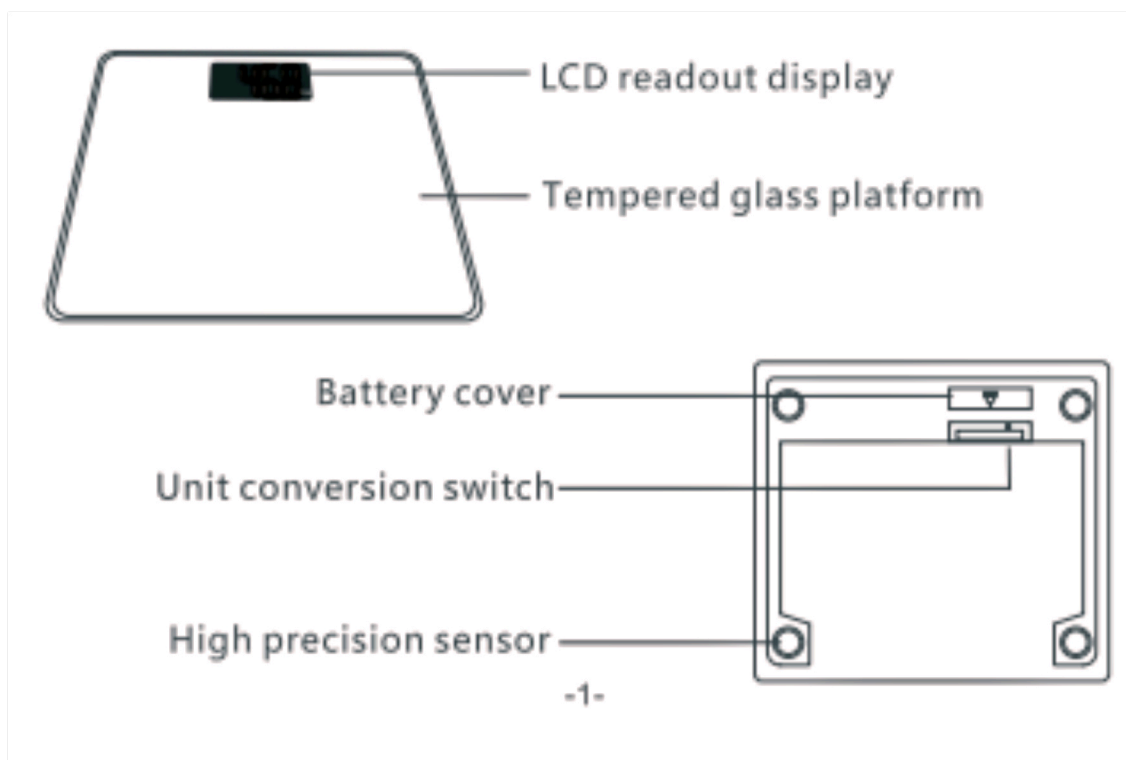Inserting or Replacing the Battery

1. On the base of the scale, gently remove the battery cover following the arrow direction.
2. Insert 2 x AAA batteries into the compartment, observing the polarity markings.
3. After inserting the batteries, replace the battery cover and slide it back into place.

Note:

Ensure the correct polarity to avoid permanent damage to the electronic components.

Operation

1. Place the scale on a solid, level surface.
2. Step gently onto the scale, and it will turn on automatically. Once a stable reading is taken, the display will lock and flash to indicate the final weight.
3. The scale will automatically turn off approximately 8 seconds after no further action.

Note:


- Do not place the scale on carpets or soft surfaces, as this will affect accuracy.
- Stand still during weighing to avoid instability affecting the results.


Unit Conversion


If the scale supports both Metric (kg) and Imperial (lb) units, you can select between them.

1. Open the battery compartment.
2. Locate the unit switch button inside. While in weighing mode, press the button to select your desired unit.

Error Indications

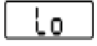 The symbol "Lo" on the display means the battery is low. Please replace it with a new battery.

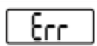 The symbol "Err" indicates that the scale is overloaded. Please step off to avoid damage.

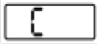 If the symbol "C" appears after stepping off, an error occurred. Step back on and weigh again.

Care and Maintenance

- Keep the scale on a flat surface, away from wet or hot conditions.
- Avoid dropping or striking the scale to prevent damage.
- Do not place heavy objects on the scale when not in use.
- Overloading the scale can damage the sensor.
- Handle the scale carefully as it is a precise electronic instrument.
- Do not immerse the scale in water.
- Clean the surface with a soft, damp cloth—avoid chemicals or harsh cleaners.
- Keep the surface dry and avoid use on slippery floors.
- Do not weigh yourself with wet feet to prevent slipping.
- This product is not for commercial or trade use.

Battery Advice

- Remove the batteries if the scale will not be used for extended periods to avoid damage from leakage.
- Keep batteries out of reach of children.
- Do not mix old and new batteries, or batteries of different types or brands, to prevent leakage or explosion.
- Do not heat, deform, or expose batteries to fire.
- Dispose of batteries according to local recycling regulations.

Note:

Electronic waste should not be disposed of with household waste. Please recycle where possible—consult your local authority or retailer for recycling advice.

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
CE DECLARATION OF CONFORMITY

The device complies with the European Union directives on safety and electromagnetic compatibility. The complete technical file of the product, which consists of the appropriate certificates as well as laboratory measurements, is located at the company's headquarters in Thermi, Thessaloniki and is available for inspection by any supervising authority.

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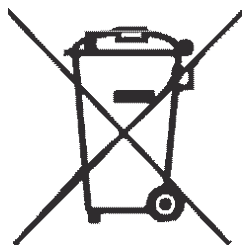Safe Appliance Disposal

This symbol means that the appliance should not be thrown away with the rest of the household waste. Look for the special recycling bins for electrical and electronic appliances. For more information, contact your local authorities or your device's retailer

**Hellenic Recycling Agency
PRODUCER REGISTRATION NUMBER: 1349**

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